



Embrace  
因倍思

# Embrace

PRE-FILLED FLUSH SYRINGES

- 1 Helps to Prevent Catheter Related Infections.
- 2 Help to Prevent Catheter Blood Reflux.
- 3 Help to Reduce the Risks of Needle Stick Injuries.
- 4 Help to Increase Clinical Staff Efficiency.




# PRE-FILLED FLUSH SYRINGES

www.emstd.cn



## Available In:

| Syringe Size | Catalogue Number |
|--------------|------------------|
| 3 mL         | catheflush3ml    |
| 5 mL         | catheflush5ml    |
| 10 mL        | catheflush10ml   |

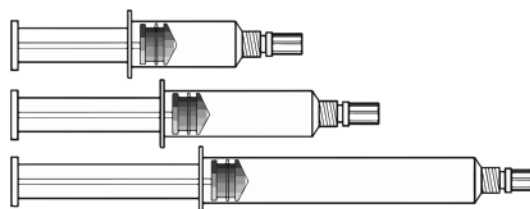

## Product Specification:

|                             |  |
|-----------------------------|--|
| <b>Product Description</b>  | Pre-Filled Normal Saline Syringe, Consist Of Core Plunger, A Syringe Barrel, A Plunger Rod, A Large Protective Tip Cap And Filled With Sterile 0.9% Sodium Chloride Injection. For <b>Single-Use</b> . It Is Suitable For Flushing All Kinds Of Vascular Access Device / Catheters.  |
| <b>Sterilization Method</b> | Steam Sterilization  |
| <b>Product Indication</b>   | <p>Embrace Normal Saline Pre-Filled Flush Syringes used for Flushing all kinds of Vascular Access Devices and Catheters to:</p> <ul style="list-style-type: none"> <li>• <b>Assess</b> catheter status/function and patency.</li> <li>• <b>Clear</b> the catheter lumen after medication administration to avoid drug incompatibilities and to ensure full dose delivery for patients.</li> <li>• <b>Lock</b> the catheters during its non-use period.</li> </ul> <p>Available in 3ml, 5ml, and 10ml sizes:</p> <ul style="list-style-type: none"> <li>• 3ml &amp; 5ml – for Peripheral Vascular Access Lines</li> <li>• 10ml- for Central Venous Catheter (CVC), Peripherally Inserted Central Venous Catheters (PICC) &amp; PortaCaths.</li> </ul> |
| <b>Package Size</b>         | 30 Each (Pre-Filled Syringe) / Box <span style="float: right;">480 Each (Pre-Filled Syringe) / Case</span>   |

## Product Features

- Specially Design Flush Device to Optimize Patient Safety.
- Individually Wrapped Syringe.
- Large Protective Syringe Tip Cap to Avoid Syringe Tip Contamination.
- Unique Double Layer Plunger Design to Reduce Blood Reflux in Catheters.
- All products - 3ml, 5ml & 10ml Syringe Sizes **Have Exact Same Standard Barrel Diameters of 10ml Syringes** to Maintain Low Pressure During Flushing Process And to Reduce The Risk of Damaging The Vein or Catheters And Comply With PICC Manufacture Flushing Recommendations.
- Clear & Bold Printed Label Graduation Scale, Also Clearly Labeled Mentioning Product Name, Syringe Content And Intended Use, To Minimizes its Misuse or Error.

• Preservatives Free

• Pyrogen Free

• DEHP Free

• BPA Free

• PVC Free

• Not Made With Natural Rubber Latex



江苏因倍思科技发展有限公司  
Jiangsu Embrace Science & Technology Development Co., Ltd.

29-2 Yingbin Road,  
Huai'an Economic and Technological Development Zone,  
Jiangsu

Version Number: ES-SA-9.0-01